

Automated Haematology Analysers

XN-9100 Series

Complex Needs,
Unlimited Solution



STANDARDISATION

EFFICIENCY

INFORMATION

UPGRADABILITY



Haematology laboratories care about the system they use because what they do is important to the welfare of the patient. Laboratory results provide insight into the clinical well-being of patients, which leads to accurate diagnosis and delivery of appropriate treatment.

Laboratories depend on Sysmex for the technology, capabilities and reliability to give them confidence.

Choosing a haematology solution is not an easy task but the Sysmex XN-Series simplifies the process and gives laboratories the freedom to define the ideal solution to meet their unique needs.

Complex Needs, Unlimited Solution

STANDARDISATION

Through compact automation

Laboratories seeking to improve processes through automation can achieve their goals in scalable phases.

Various functional modules can be integrated with the XN-9100 allowing automation to be introduced gradually to the laboratory. These components include:

- XN-Series automated haematology analysers with a comprehensive test menu to accommodate a variety of laboratory settings and test volumes.
- Fully automated slide maker and stainer, SP-50, that optimises the smearing conditions for every sample based on haematocrit.
- Digital morphology analyser, DI-60, captures cells and performs pre-classification automatically.
- Tube sorter, TS-10, sorts and archives samples into pre-defined sorting areas.

1

TS-10 - integrated tube sorter and archiver

Benefits:

- Allows sorting of abnormal samples for quick identification and timely action leading to a reduction in turnaround time (TAT).
- With an upgrade kit, sorting speed and capacity can be expanded. This allows laboratories to keep the initial investment and have a scalable solution to meet an increase in workload or sorting capacity.

2

DI-60 - integrated digital morphology

Benefits:

- Standardisation of differential results, operational and validation processes, delivering consistency.
- Facilitates collaboration, review and consultation between morphology experts, creating flexible and cost-effective utilisation of staff and resources through the Sysmex remote review software.
- Monitors and promotes cell morphology competence with a web-based proficiency application.
- Reduces the neck, wrist and eye fatigue associated with traditional microscopy.

3

SP-50- fully automated slide maker and stainer

Benefits:

- Standardisation and efficiency in smear preparation for timely morphology assessment and interpretation.
- Slides are only prepared when samples fulfill user pre-defined conditions, promoting workflow efficiency based on clinical needs.
- Slim design makes it ideal for labs with limited space.

1

2

3

4

4

XN-Series analysers are comprised of XN-10 (green panel)

- CBC+DIFF+NRBC as standard.

XN-20 (blue panel)

- CBC+DIFF+NRBC+RET+WPC as standard.

Benefits:

- Enhanced diagnostic capability.
- Flexibility with optimisation of an integrated solution.
- Faster test results supporting clinical decisions for better patient management.



EFFICIENCY

Truly walkaway system

Laboratory technologists today face growing demands to multi-task. Expansion of duties such as sample processing, analysis of results, sample archiving and documentation can impact KPI achievements.

The XN-9100 series is a state-of-the-art scalable solution that helps address these growing demands and challenges.

With the twin XN analysis module, a leaner automated slide staining /making unit, SP-50 and on-board decision rules, this compact integration allows unparalleled efficiency and productivity for even space-constrained laboratories.

Using Sysmex technology to maximise efficiency, operators only need to load and unload samples on the XN-9100, enabling staff to focus on higher value.



Productivity comes in many different forms. With the modular design, flexible configuration and overall footprint reduction of the XN-9100 systemisation solution, the limitation lies only with you!

Enhanced efficiency and productivity through streamlined workflow and onboard rule set

With onboard rerun and reflex capability, the XN-Series can automatically reanalyse samples with abnormal or unreliable results as per the defined laboratory criteria. The benefits include:


- Delivery of highly reproducible results in the shortest time possible improving turnaround time.
- Minimal manual interventions for a more efficient use of time and manpower resources.
- Standardised sample processing workflow, promoting peace of mind in the laboratory.

Enhanced automation workflow through flexible systemisation layout

The XN-9100 Series comes in five alternative layouts to accommodate existing building fixtures such as pillars or water supply points without the laboratory undergoing any major renovations.



Single master switch for easy startup of the XN-9100

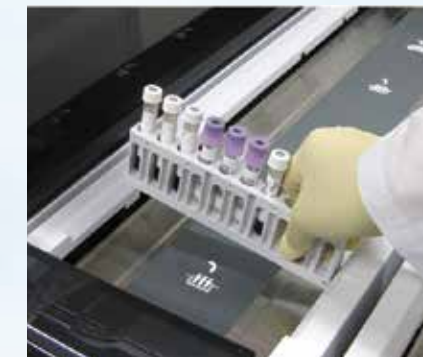

Continuous analysis from sample loading to reflex testing to sample archiving



Automatic sample analysis, including reflex testing, according to LIS test orders



Enhanced workflow with an integrated smear preparation and cell-pre-classification unit based on user-defined rules and conditions



INFORMATIVE

Trustworthy technology providing insight on patient well-being

The XN-Series analysers utilise laser flow cytometry to enumerate blood cells. The complex algorithms are applied to the measured cellular characteristics for classification of WBCs, RBCs, PLTs as well as for flagging of abnormal cellular populations.

Numerous parameters are also available providing information about different disease states:

1 Response to erythropoietic stress

NRBC (Nucleated Red Blood Cells) are a useful indicator of erythropoietic stress or red cell related disorders such as thalassemia, myeloproliferative disorders etc. The duration and presence of NRBCs is also associated with poor prognosis in critically ill patients¹.

On XN-Series

- NRBC is a standard parameter with every CBC analysis.
- Stained directly with a linearity up to 600 NRBC/100 WBCs.
- Automatic correction of WBC when NRBCs are present.

Benefits:

- Increase productivity and improve sample TAT.
- Cost reduction as additional testing for NRBCs is no longer necessary.
- Therapeutic support with reliable NRBC counts.

2 Supports infection/ inflammation diagnosis and therapy monitoring

Immature Granulocytes (IG) supports differentiation between SIRS and sepsis and a useful marker in monitoring infection/ inflammation in patients undergoing therapy².

On XN-Series

- IG is a standard parameter with every DIFF analysis, providing a 6-part differential.
- IG includes metamyelocytes, myelocytes and promyelocytes.

Benefits:

- Enhance workload efficiency through reduction in blood smear and manual microscopy.
- Support diagnostic and prognostic use in infection/inflammation and therapeutic monitoring in combination with other parameters such as cytokines.

3 Differentiates hypo-production from peripheral destruction as a cause of thrombocytopenia

Immature Platelet Fraction (IPF) are young, reticulated platelets containing residual RNA. IPF reflects thrombopoiesis with a low count pointing to impaired platelet production, while a high IPF count suggests peripheral destruction. It is an early indicator of bone marrow recovery and a potential marker for a more optimised platelet transfusion approach^{3,4}.

On XN-Series

- Fluorescent platelets and IPF are both reported from the PLT-F channel where platelet-specific stain is used and cells are counted 5x more.
- Accurately reflects the rate of thrombopoiesis.

Benefits:

- Aids in appropriate treatment decision.
- Supports optimised platelet transfusion approach.

4 Accurate timing of peripheral blood stem cell transplant (PBSCT) harvest

High comparability between **Human Progenitor Cells (HPC)** measurement and CD34 analysis supports application in PB stem cell collection.⁵

On XN-Series

- Human progenitor cells (HPC) are counted in a dedicated mode of measurement in the WPC channel.
- Automated analysis with good comparability with CD34+cells.

Benefits:

- Rapid analysis allows for a prompt determination of optimal stem cell harvest time.
- Realise savings in workflow and resources.

5 Body fluid (BF) analysis in infection or inflammatory conditions

The Body Fluid (BF) analysis in the XN system provides rapid and accurate measurement of body fluids such as the CSF, synovial fluid, pleural fluid and CAPD. The 2-part differential results of mononuclear (MN) and polymorphonuclear (PMN) cell population aid in the quick distinction between viral and bacterial infection.⁶

On XN-Series

- Automated body fluids analysis is done in a dedicated mode of measurement.
- No sample pretreatment and additional reagents are required.

Benefits:

- Improve sample turnaround time (TAT)
- Better standardisation and reproducibility of results reducing errors commonly associated with manual procedures.
- Quick distinction between viral and bacterial infection guides clinical decisions for patient treatment.

6 Diagnosis and assessment of iron status and therapy.

Ret-He (Reticulocyte Haemoglobin Equivalent) The haemoglobin content of reticulocyte, reflects real-time iron bioavailability and quality of erythropoiesis. A low Ret-He value means iron is lacking or is not bioavailable for erythropoiesis. The Ret-He parameter differentiates the classical cause of iron deficiency from a functional cause. It can also be used to monitor response to EPO or iron therapy.⁷

On XN-Series

- Ret-He is available with every retic analysis along with immature retic fraction (IRF), Hypo-He and Hyper-He.
- Hypo-He and Hyper-He together with %Micro,%Macro, are also useful in classification of various types of anaemia.

Benefits:

- RET-He helps differentiate between functional and classical iron deficiency,
- Ret-He, in combination with other retics-related parameters, allows clinicians to draw conclusions on the quality and quantity of the young RBC fraction supporting effective treatment decisions.
- Ret-He aids in monitoring responses to EPO and IV iron therapy.

UPGRADABILITY

Modularity concept that caters to changing needs

Flexible activation of applications



Applications

CBC+DIFF



NRBC increases in response to erythropoietic stress.
 IG aids in the diagnosis, prognosis and therapy monitoring of infection/inflammation.



+ Body Fluids application



Automated two-part differential body fluids results allowing quick distinction between viral and bacterial infection.



+ RETICS & PLT-F applications



RET-He and IPF provides insight about the causes of anaemia and thrombocytopenia.



+ HPC application



HPC helps to assess optimal timing of peripheral blood stem cell transplant (PBSCT) harvest.

The XN-Series at a glance

The Standard Diagnostic Components for each XN Module

Analytic Components: XN-CBC and XN-DIFF are standard.
All other applications are optional




XN-1000
First step into full automation



Hourly throughput:
Up to 100 samples an hour

- Onboard decision rules with user-defined rerun/reflex capabilities
- Customisable clinical applications to cater to variable clinical needs

XN-2000
Workload optimisation

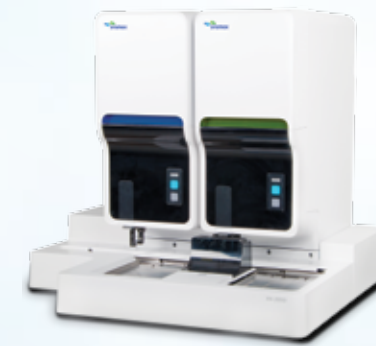

Hourly throughput:
200 samples an hour

- Unique co-primary solution
- Automatic workload balancing between the two analysers
- Reagent sharing option is available

XN-1500
Simplifying Automation



Hourly throughput:
100 CBC samples with 30-75 slides an hour

- Automates sample analysis and slide preparation
- Reduces manual intervention through auto-rerun/reflex capability
- Compact solution suitable for most laboratories

XN-9100
Scalable Automation, Compact Solution

- Broad capabilities catering to different clinical demands
- Workflow optimising configurations
- Flexible, modular and compact integration



XN-3100
Compact Automation



Hourly throughput:
200 CBC samples with 30-75 slides an hour

- Scales up automation for medium-sized laboratories
- Standardised and seamless workflow from CBC analysis to integrated slidemaking/staining unit
- Enhances efficiency through onboard reflex rule set
- DI-60, Digital Morphology integration is an option

NOTES

Disclaimer

The uses or clinical applications described are based upon published scientific and clinical evidences.

It is the clinician's responsibility to validate any off-label applications for use in routine clinical practice.

Notice of Intended Use

The XN-Series modules are quantitative multi-parameter automated hematology analysers intended to be used in clinical laboratories for *in vitro* diagnostic use in screening patient populations.

The publication of this material is sponsored by Sysmex Asia Pacific for health care professionals. Reproduction of this material in whole or part is not permitted without the permission of Sysmex Asia Pacific.

References

1. Monteiro Júnior JGdM, Torres DdOC, da Silva MCFC, Ramos TMdB, Alves ML, Filho WJN, et al. (2015) Nucleated Red Blood Cells as Predictors of All-Cause Mortality in Cardiac Intensive Care Unit Patients: A Prospective Cohort Study. PLoS ONE 10 (12): e0144259.
2. Nierhaus et al.: Revisiting the white blood cell count: immature granulocytes count as a diagnostic marker to discriminate between SIRS and sepsis - a prospective, observational study. BMC Immunology 2013 14:8.
3. Morkis IVC, Farias MG, Rigoni LDC, Scotti L, Gregianin LJ, Daudt LE, et al. Assessment of immature platelet fraction and immature reticulocyte fraction as predictors of engraftment after hematopoietic stem cell transplantation. Int J Lab Hematol. 2015; 37, 259–264
4. Briggs C, Hart D, Kunka S, Oguni S, Machin SJ. Immature platelet fraction measurement: A future guide to platelet transfusion requirement after haematopoietic stem cell transplantation Transfus Med. 2006; 16, 101–109
5. Peerschke EI, Mounq C, Pessin MS, Maslak P. Evaluation of new automated hematopoietic progenitor cell analysis in the clinical management of peripheral blood stem cell collections. Transfusion. 2015 55(8): 2001–2009.
6. Fleming CKA. From manual microscopy to automated cell counters for first line screening of body fluids: "But not without a special body fluid mode." Ned Tijdschr voor Klin Chemie en Lab. 2016; 41: 229-234
7. Thomas C, Kirschbaum A, Boehm D, Thomas L. The diagnostic plot: A concept for identifying different states of iron deficiency and monitoring the response to epoetin therapy. Medical Oncology. 2006 (23), 1:23–36

Registered under Act 737
Registration No: IVDB17531114618

Sysmex Asia Pacific Pte Ltd
Tel +65 6221-3629 Fax +65 6221-3687
www.sysmex-ap.com

Sysmex India Pvt. Ltd
Tel +91 (22) 6112-6666
www.sysmex.co.in

PT Sysmex Indonesia
Tel +62 (21) 3002-6688 Fax +62 (21) 3002-6699
www.sysmex.co.id

Sysmex (Malaysia) Sdn Bhd
Tel +60 (3) 5637-1788 Fax +60 (3) 5637-1688
www.sysmex.com.my

Sysmex New Zealand Ltd
Tel +64 (9) 630-3554
www.sysmex.co.nz

Sysmex (Thailand) Co., Ltd
Tel +66 (2) 032-2536 Fax +66 (2) 116-5396
www.sysmex.co.th

Sysmex Philippines Inc.
Tel +63 (2) 621-2460 Fax +63 (2) 621-2432
www.sysmex.com.ph

Sysmex Vietnam Co., Ltd
Tel +84 (8) 3997-9400 Fax +84 (8) 3997-9405
www.sysmex.com.vn

Sysmex Vietnam Co., Ltd (Hanoi Branch)
Tel +84 (4) 3776-7020 Fax +84 (4) 3776-7022
www.sysmex.com.vn

Sysmex Australia Pty Ltd
Tel +61 (2) 9016-3040
www.sysmex.com.au

Design and specifications may be subject to change due to further product development.